

Pressure / Vacuum Meter

DigiMano 1000



Features:

- Accurate
- Pocket Size
- Dependable
- Luer Lock Fitting
- Large LCD Display
- Two Engineering Units
- Low Battery Indication
- Liquid and Gas Compatible
- 9 Volt Battery Operated
- Measures Pressure and Vacuum

The **DigiMano 1000**, Pressure Vacuum Meter is a pocket size digital instrument designed to measure the pressure and vacuum of non-corrosive fluids and gases. The versatile instrument displays the reading on a large 3½ digit LCD Display

The **DigiMano 1000**, Pressure Vacuum Meter has an accuracy of 0.25% full scale. Seven models with different engineering units and pressure ranges are available.

To provide superior performance, the **DigiMano 1000**, Pressure Vacuum Meter incorporates one of the latest semiconductor sensors. It is powered by one 9V alkaline battery or an AC Adapter (Optional)

DigiMano 1000 Shown with Syringe Pump



Specifications:

Display: Large LCD 3 ½ digits

Measures: Pressure or Vacuum.

Accuracy: 0.25% Full Scale \pm 1 Digit

Over Pressure: Three Times Pressure Range

Medium Compatibility: Non-Corrosive fluid and gas


Eng. Unts : A combination of units are available depending on the Model.

Models / Ranges: Several Models to Choose from ranges 2PSI to 20 PSI

Power: 9 Volt Alkaline Battery or Optional AC Adapter

Dimensions: 3 x 2 x 1.4 inches (7.6 x 5 x 3.6 cm)

Weight: 0.2 lbs **Analog Output:** 3.5mm Sterio Jack (0-2 Volt)

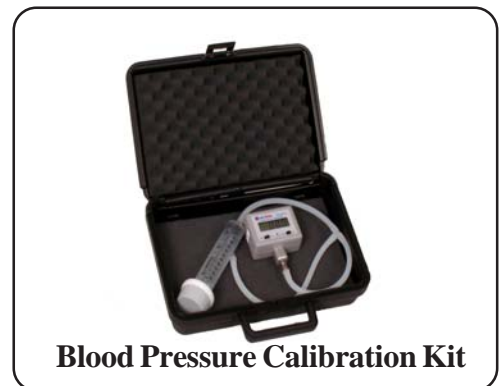

Blood Pressure Calibration Kit

Ordering Information:

Part #	Range / Units	Part Number	Range / Units
200-200 IN:	199.9 mmHg/199.9 inH2O	200-2000 IN:	1999 mmHg/199.9 inH2O
200-200 PS:	199.9 mmHg/19.99 PSI	200-20 PS:	1999 mmHg/19.99 PSI
200-200 CmH20:	199.9 cmH20/19.99 PSI	200-200 kPa-InH20:	199.9 kPa/199.9InH20
200-200 kPa:	199.9 kPa/199.9 cmH20	204	Syringe Pump



110 Toledo Street, Farmingdale, NY 11735, Phone: 631-531-0100, Fax: 631-531-0101

www.Netech.org 1-800-547-6557